


Blister Pack Opener

USER QUICK GUIDE

Introduction

The Blister Pack Opener is designed to help users with arthritis or limited hand strength open blister packs of pills.

Features

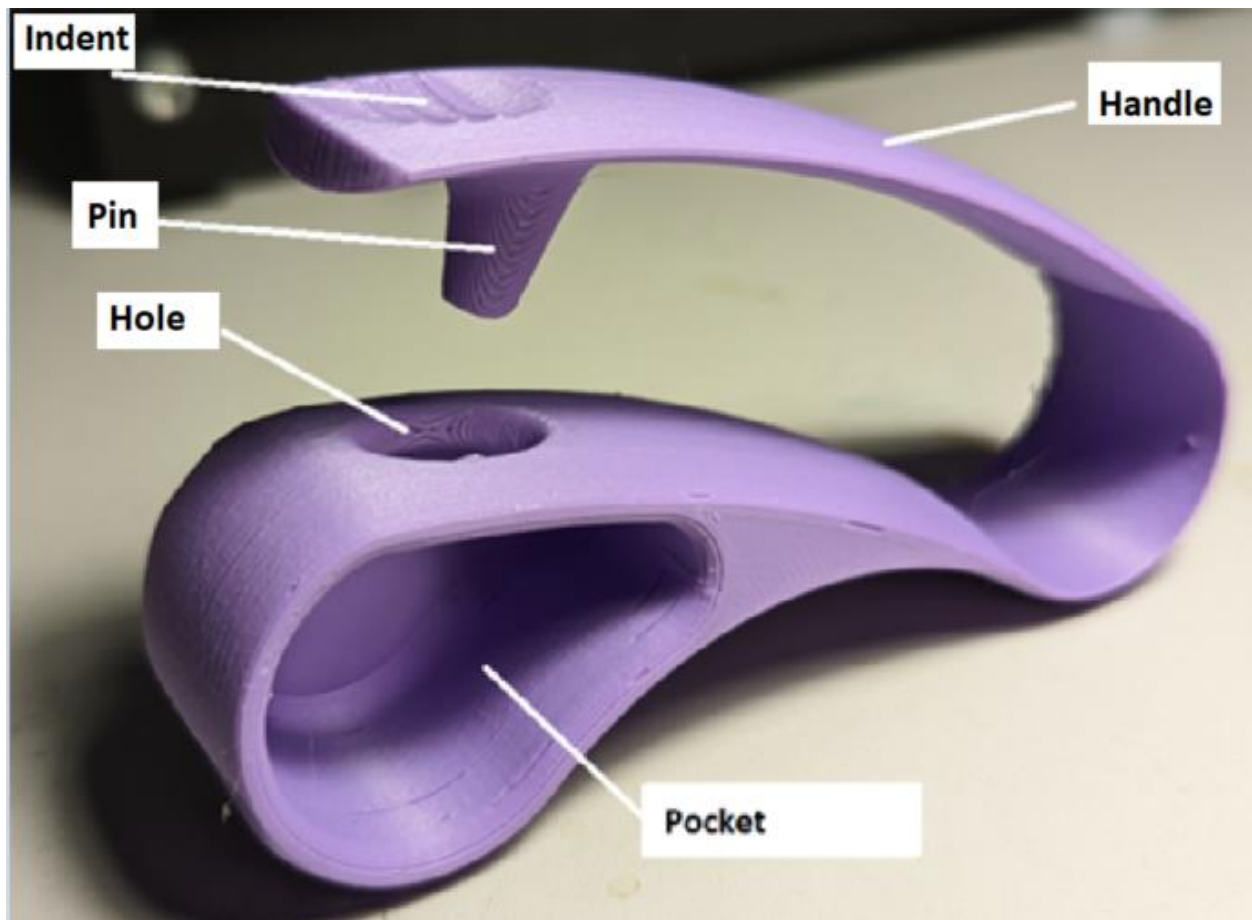

Usage

To use, place the blister pack so the flat side is over the hole, and the pill is facing the pin. Squeeze the handle so that the pill is popped out of the blister pack so that the pill pops out into the tube in the handle.

Blister Pack Opener

USER QUICK GUIDE



Compatibility

The device works with any pill that fits in a 20x15mm oval.

Specifications

Length: 10cm

Width: 6cm

Height: 2.5cm

Weight: 20g

Cleaning

The device can be wiped clean with a damp cloth.